



Technology deployments

Boost your capabilities with new technology deployments

EasiPak technology deployments provide the right software and clinical training to enhance your team's capabilities











Scanner connectivity

An interface that allows communication between your scanners and the MEDRAD® Stellant CT Injection System



Enhances clinical capabilities by synchronizing scan timing and simplifying workflows.

- Provides precision scanning techniques and improves contrast bolus timing
- Can save time by starting the procedure from the control room
- Limits the amount of back and forth attention between injector and scanner controls
- Allows the flexibility to remain at your patient's side









Scanner connectivity

Use one button to simultaneously trigger the scanner and the injector




New process



Technologist stays by patient's side and monitors patient during initial injection

One button on either the scanner or the injector can initiate both the injection and scanner protocol



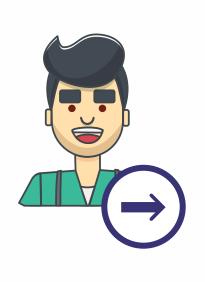

Technologist sets up patient and injector protocol


Technologist walks to scan room

> Technologist sets up scanner protocol

> Technologist walks back to control room





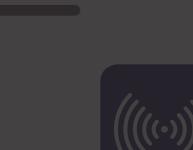

Scanner connectivity



New process



Technologist stays by patient's side and monitors patient during initial injection

One button on either the scanner or the injector can initiate both the injection and scanner protocol

Current process



Technologist sets up patient and injector protocol

Technologist walks to scan room

Technologist sets up scanner protocol

Technologist walks back to control room

Scanner connectivity is available for select SIEMENS, TOSHIBA, GE, PHILIPS and HITACHI scanners.





BAYER E R

Software designed to automatically capture and distribute key data for contrast enhanced procedures



See how Automated documentation reduces administrative tasks

Allows radiology teams to reduce administrative rework and multiply efficiencies.

- Import patient and procedure details from Modality Worklist
- (Aa) Auto-populate individual contrast injection records into PACS and Speech Recognition
- Improve reporting accuracy and traceability
- Reduce dictation and rework time for Radiologists








Select patient directly from Modality Worklist/Patient Worklist

Patient Worklist						
	Patient Name	DOB	ID	Time & Date	Accession Number	Study Description
<u>i</u>	SEXTON, EMMANUEL	05/06/2007	ID-00000	13:48 PM 09/22/2011	SSS-000000000806	Study info goes here: 0806
i	PECK, RUFUS	12/01/1990	ID-00000	14:48 PM 09/22/2011	SSS-000000000061	Study info goes here: 0061
(i)	ATKINS, LINO	01/05/1992	ID-00000	14:48 PM 09/22/2011	SSS-000000000291	Study info goes here: 0291
(i)	MORSE, ALPHONSE	09/05/1976	ID-00000	14:48 PM 09/22/2011	SSS-000000000726	Study info goes here: 0726
(i)	HARPER, RALEIGH	08/27/1993	ID-00000	15:48 PM 09/22/2011	SSS-000000000164	Study info goes here: 0164
i	MIDDLETON, BARTON	02/08/1993	ID-00000	15:48 PM 09/22/2011	SSS-000000000984	Study info goes here: 0984
i	MEJIA, STEPHAN	11/15/1955	ID-00000	16:48 PM 09/22/2011	SSS-000000000055	Study info goes here: 0055
i	MUNOZ, EMILE	08/19/1964	ID-00000	16:48 PM 09/22/2011	SSS-00000000180	Study info goes here: 0180
i	VANG, DALLAS	10/29/1920	ID-00000	16:48 PM 09/22/2011	SSS-000000000478	Study info goes here: 0478
i	SUMMERS, EDDIE	06/02/1979	ID-00000	16:48 PM 09/22/2011	SSS-000000000529	Study info goes here: 0529
Patient Name: Update Cancel						

All patient data that appear in this document are fictitious. No actual patient information is shown.



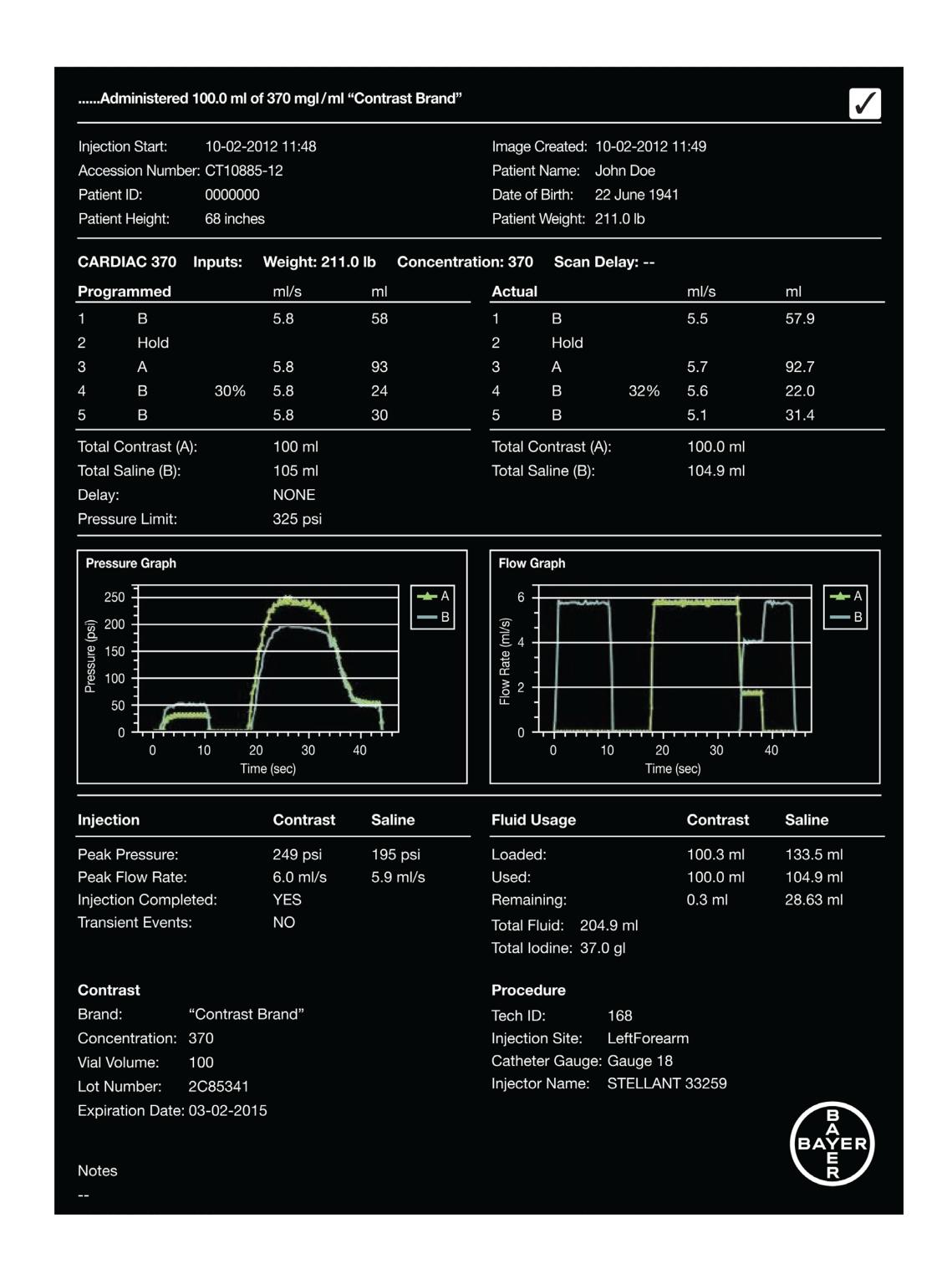


Enter and capture vital injection details directly into PACS

All patient data that appear in this document are fictitious. No actual patient information is shown.



Enter and capture vital injection details directly into Speech Recognition

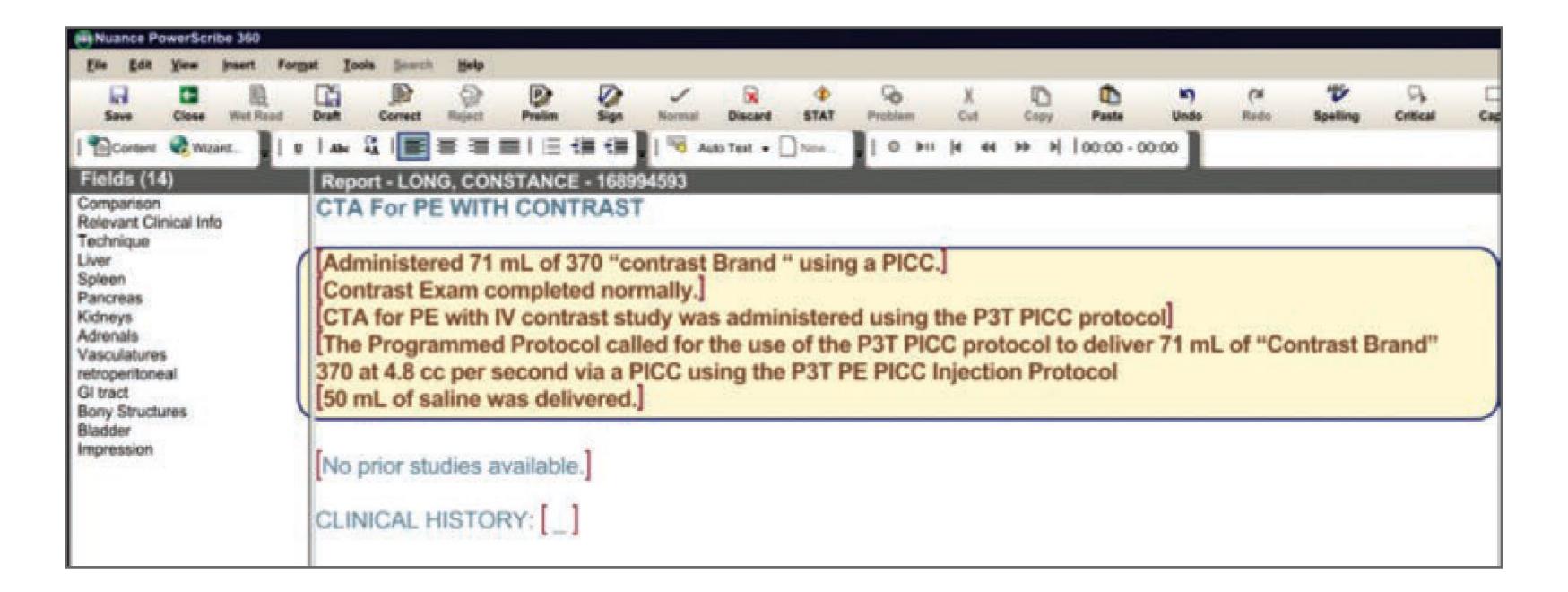


All patient data that appear in this document are fictitious. No actual patient information is shown.





Complete the same data capture and reporting process in fewer steps

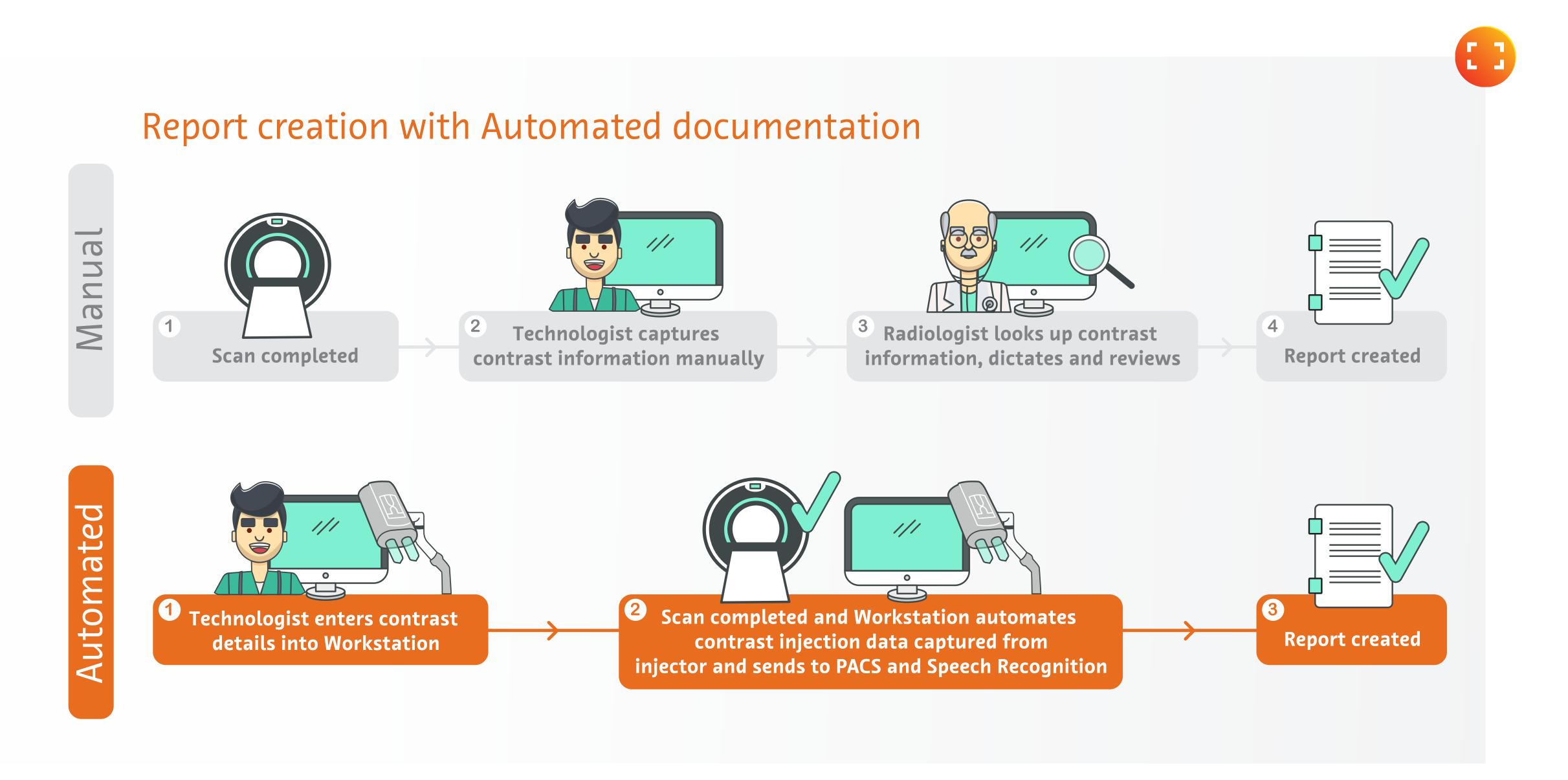






BAYER E R

Automated documentation

Complete the same data capture and reporting process



Report creation with Automated documentation

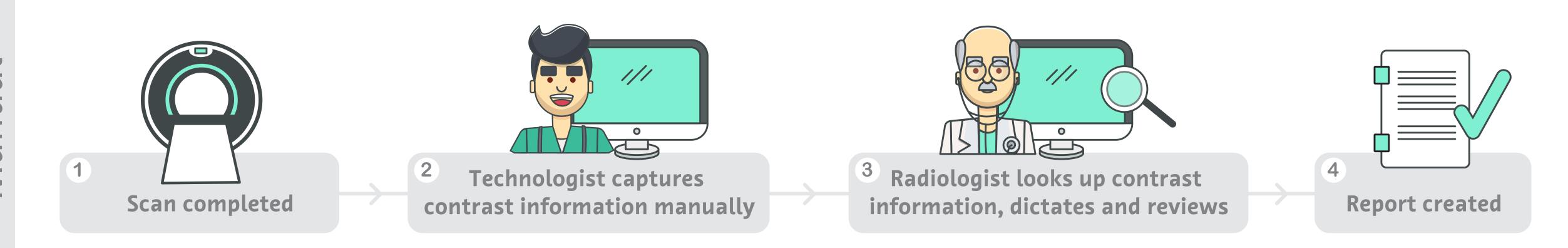














Personalized protocols

An algorithm that allows technologists to personalize injection protocols and provide consistent image quality amongst clinicians

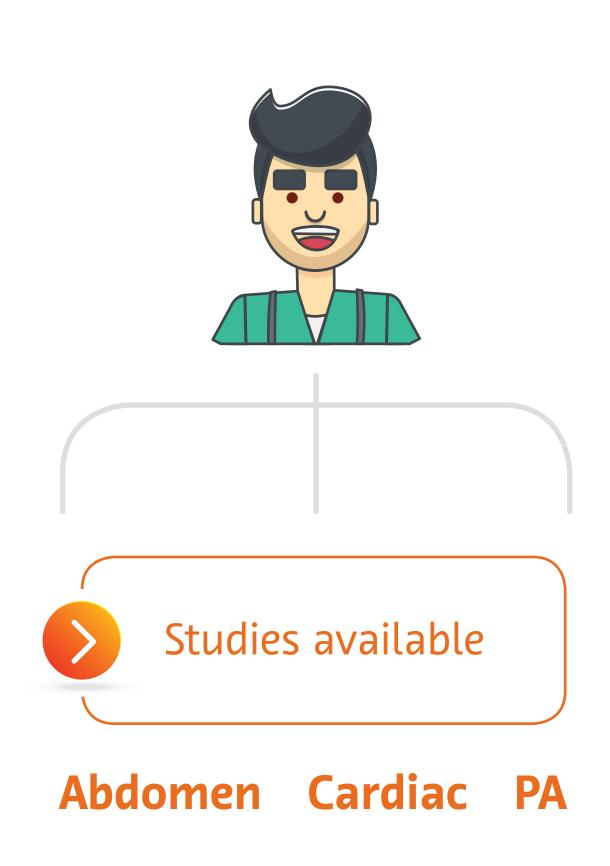

Allows technologists to account for unique patient and study variables.

- Calculate contrast media protocols automatically
- Personalize the volume and flow rate for each patient
- © Configure according to the workflow, policies and preferences of each facility








Personalized protocols


A clinically validated algorithm designed to provide consistent image quality amongst clinicians by managing personalized contrast delivery protocols with flexibility to incorporate pre-set Radiologist preferences for dosing, contrast concentrations and scan timing

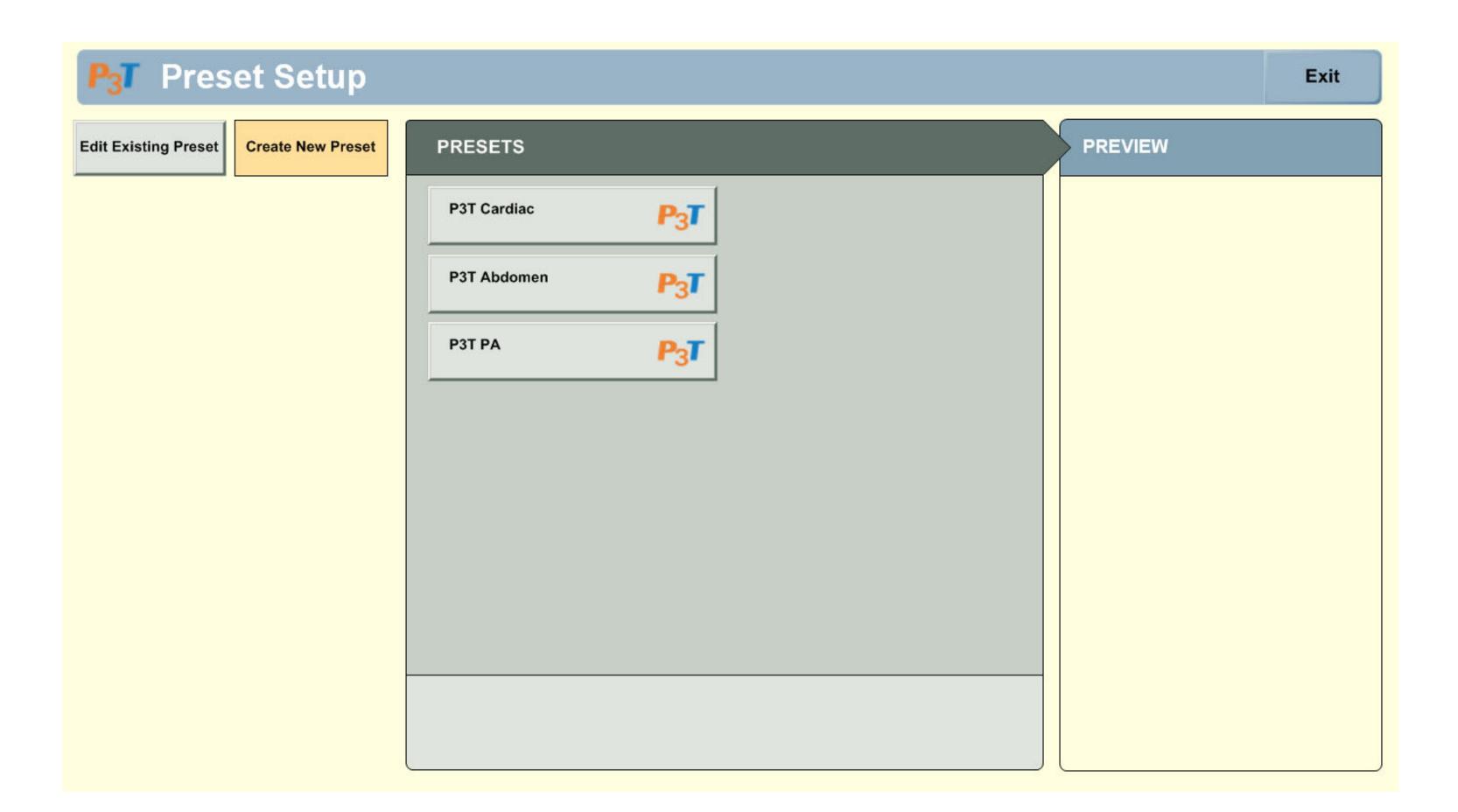



Cardiac studies


Manage cardiac-imaging complexity

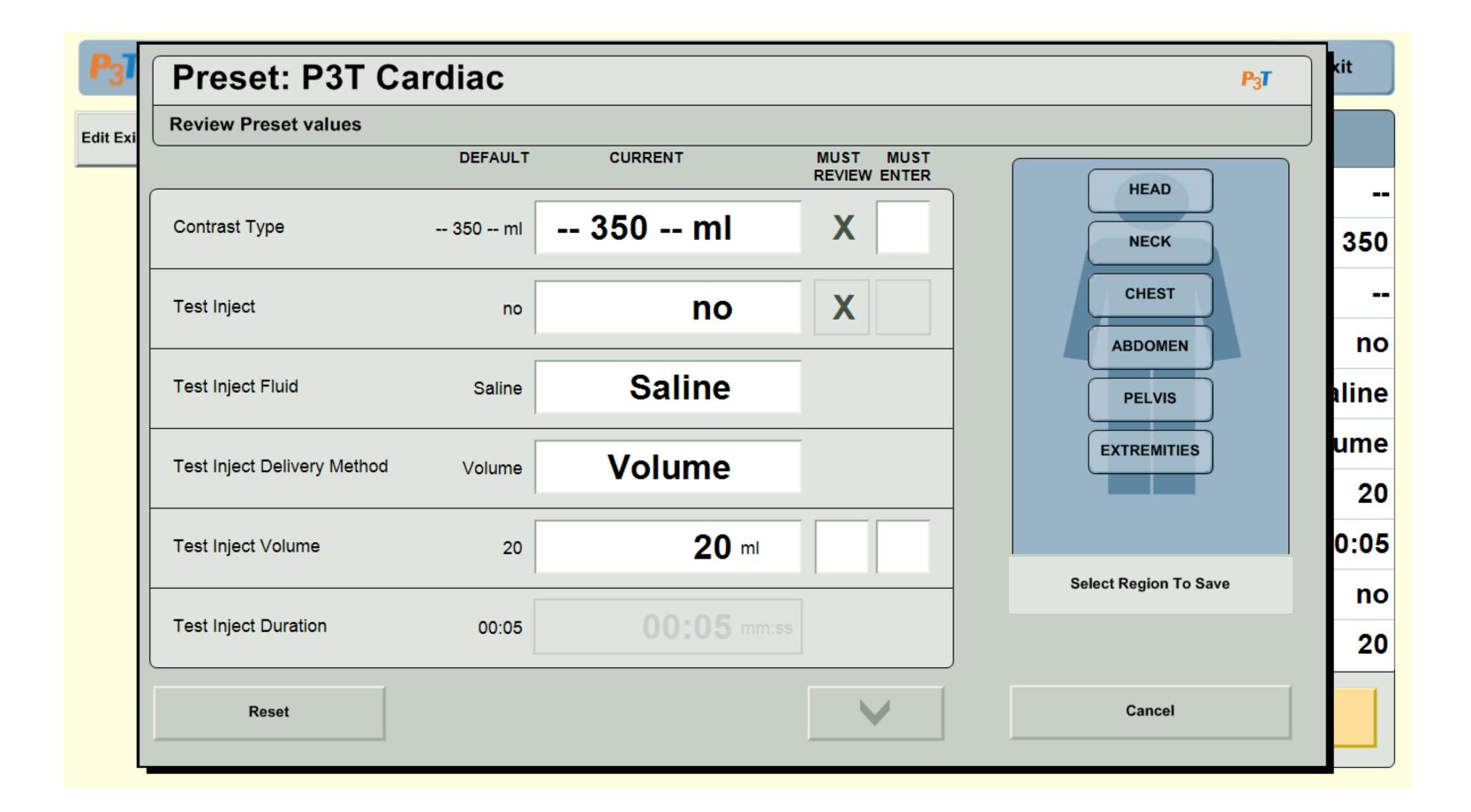



Abdomen studies


Facilitate implementation of multiple CT abdomen protocols

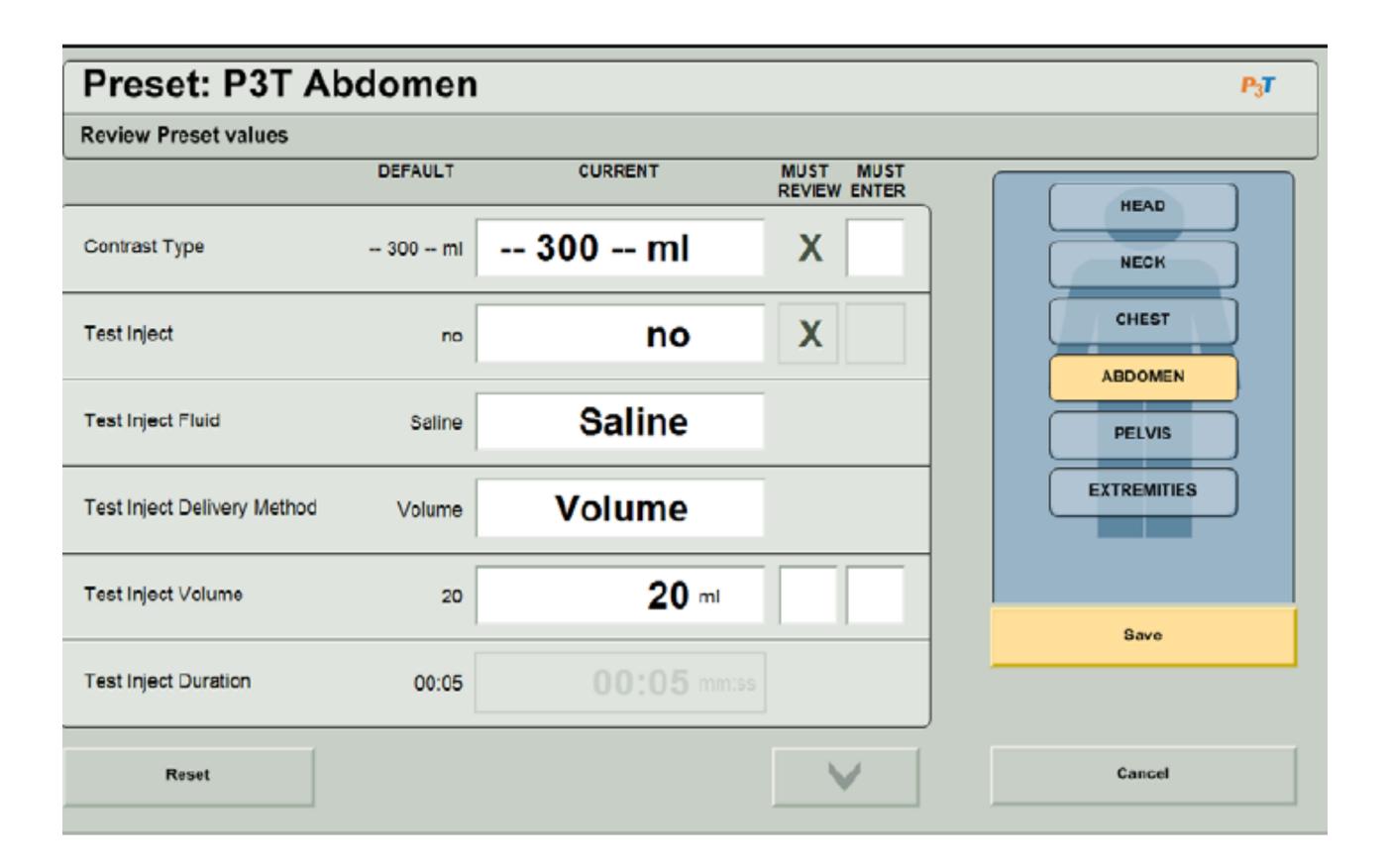



Pulmonary angiography studies

Reducing occurrences of suboptimal CT pulmonary angiography studies, albeit at a slightly higher contrast dose than the standard protocol of 80 mL¹

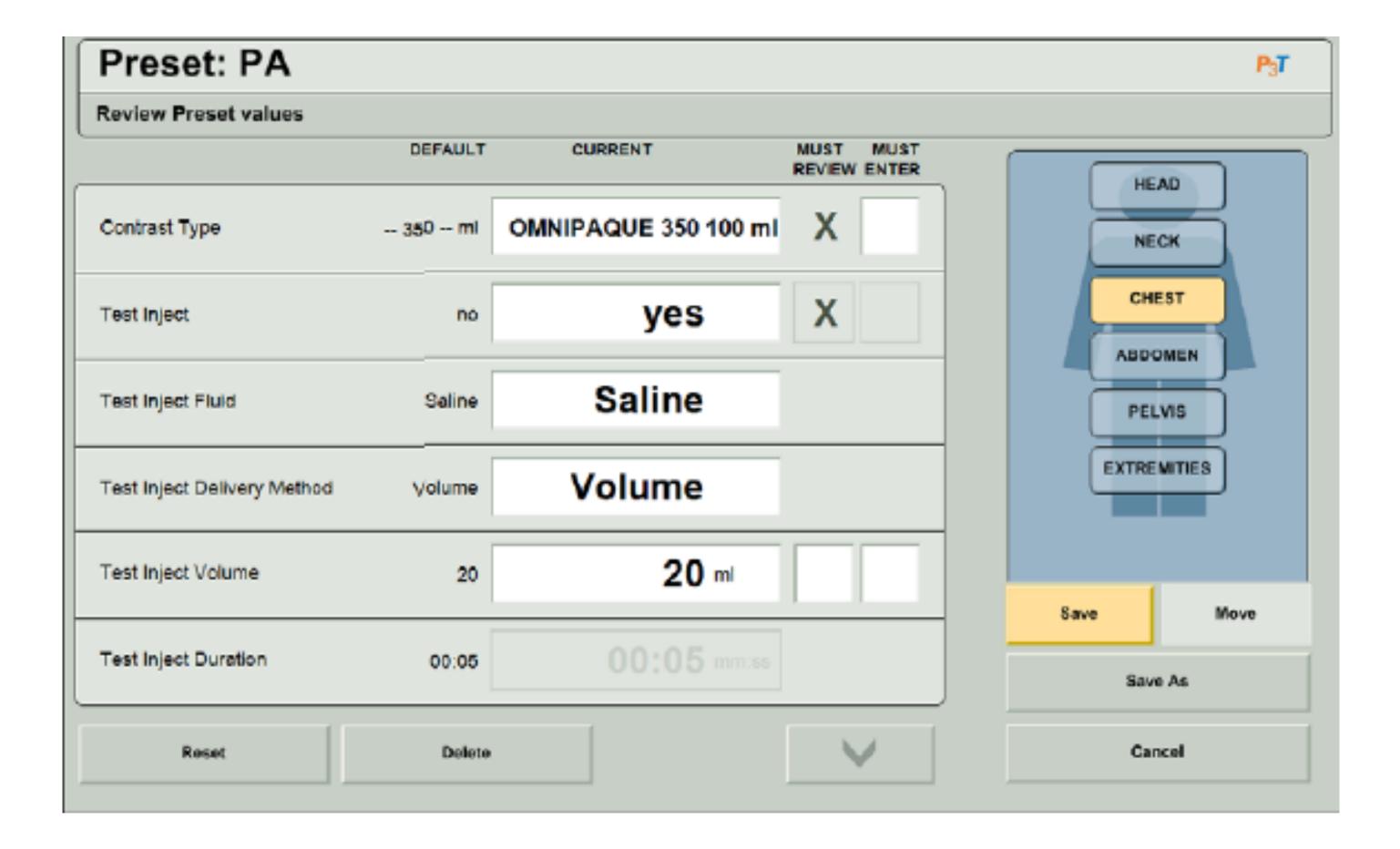


Reference: 1. Deible CR, et al. A Clinical Evaluation of an Automated Software Program (Certegra® P3T® PA) for Patient Specific Contrast Injection During Chest CTA to Exclude Pulmonary Embolism. RAD-INF-14-06087 May 2014.



Personalized protocols

Edit an existing pre-set

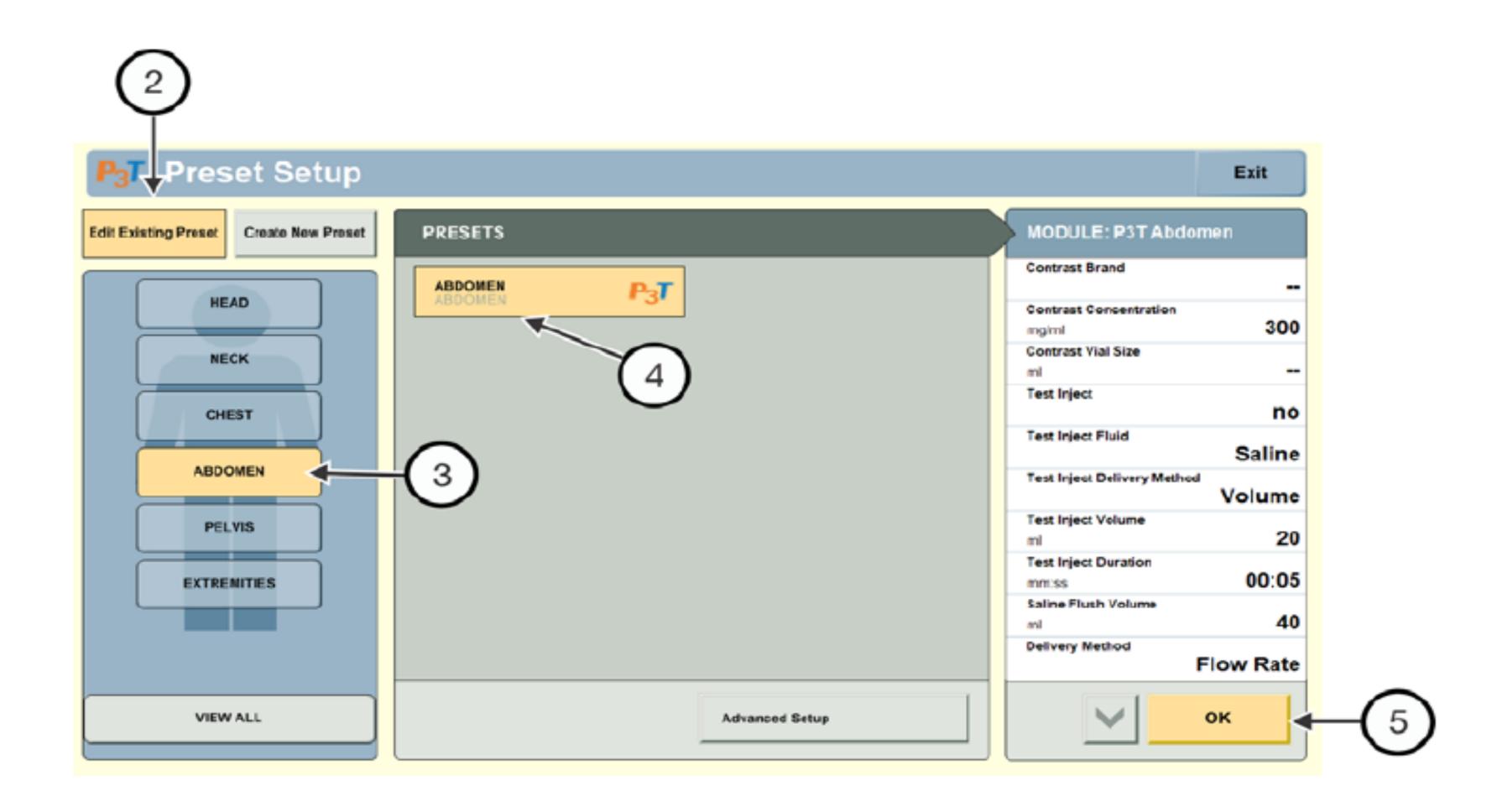






Data-driven insights

A platform that connects contrast, injector and scan information to provide critical insights

Allows radiology practices to generate diagnostic quality images more safely, consistently and efficiently.

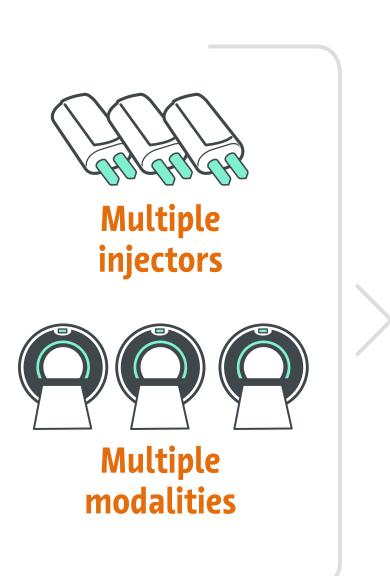








quality

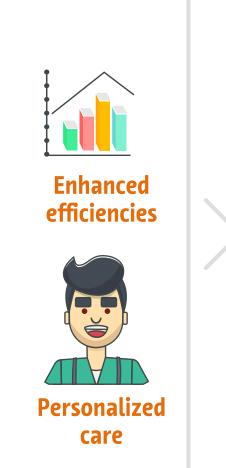

- Easily track injection volumes, brand and vial size by injection protocol, device or across the enterprise
- Track injection protocol adherence to optimize appropriate injections
- Compare point-of-care data capture adherence across the enterprise for quality management
- Analyze scanner throughput for contrast enhanced exams

Future technology deployment



